

MedLinket


Speciality for getting and transferring signals!



AM-601 Pulse Oximeter



Patent Algorithm

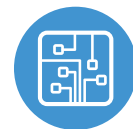

Imported Chip



OLED Screen



8s Automatic Power Off




Screen Switch

- The alarm upper and lower limits of blood oxygen saturation and pulse rate can be set, and the alarm will be given automatically if the limit is exceeded.
- Switchable between displays, optional waveform interface and large character interface
- Patented algorithm, accurate measurement under weak perfusion and jitter
- OLED display, no matter day or night, can display clearly
- Low power, long battery life, low use cost

* The picture is for reference only, the actual item is the standard.

Pulse Oximeter

Listen To Fingertips Heartbeat, Feel Temperature of Life

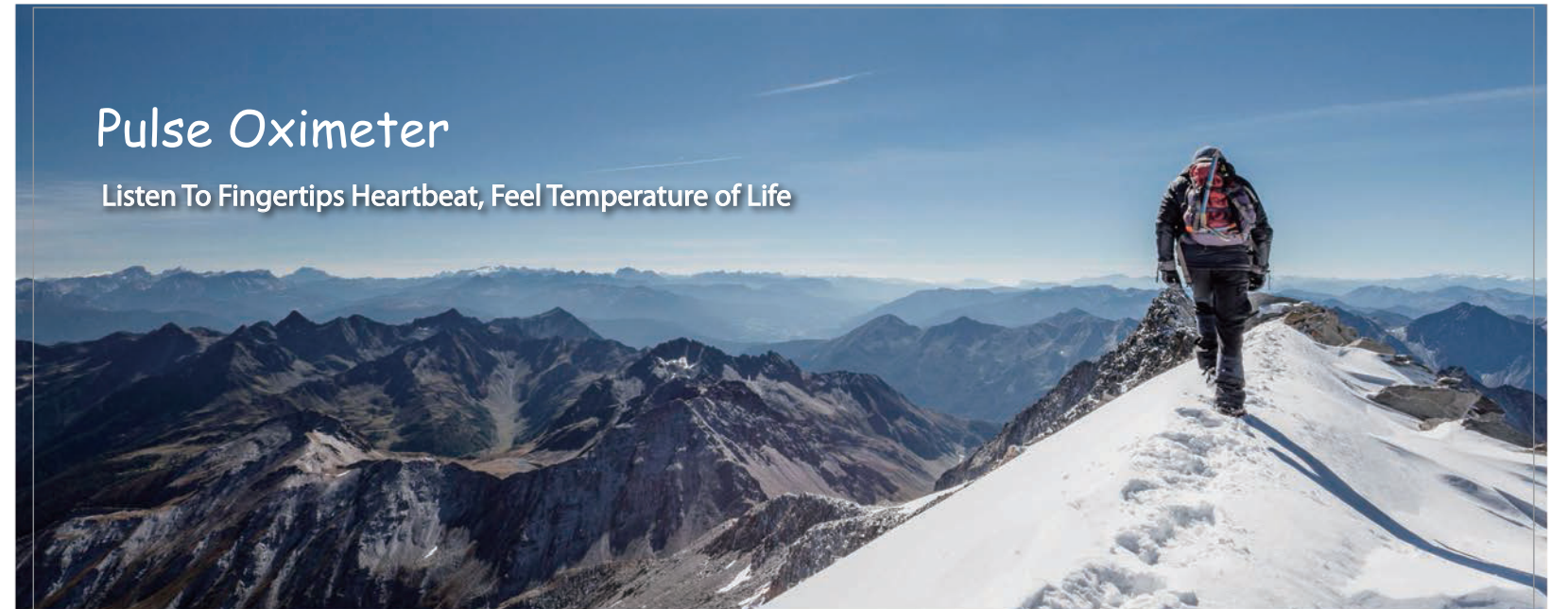

Basic Specification

Product Name	Pulse Oximeter	Product Model	AM-601
Display Screen	OLED Screen	Display Direction Switch	Display 4 Directions Switch
Weight/Size	31.5g/L*W*H : 61*34*30.5 (mm)	Automatic Alarm	Available For Setting Upper & Lower Limit, Automatic Alarm When Beyond Limit
Measuring Range	SpO2 : 70 ~ 100% Pulse Rate : 30 ~ 250bpm	Measuring Display Unit	SpO2:1% Pulse Rate:1bpm
Power	DC 3.0V (2 Pieces AAA alkaline batteries)	Measuring Accuracy	SpO2 : 90% ~ 100% , ±2% ; : 70% ~ 89% , ±3% ; Pulse Rate : ±3bpm
Accessories	English user manual / Lanyard	LED Wave Length	Red Light : About 660nm ; Infrared Light : About 905nm

Shenzhen Med-link Electronics Tech Co., Ltd

4th and 5th Floor, Building Two, Hualian Industrial Zone, Xinshi Community, Dalang Street, Longhua District, 518109 Shenzhen, China

Fax: +86 755 61120055 E-mail: kefu@medxing.com

Website: www.medxing.com



Med-linket
Official Website